



Tera Harz Care

ULTRASONIC CLEANER AND WARMER

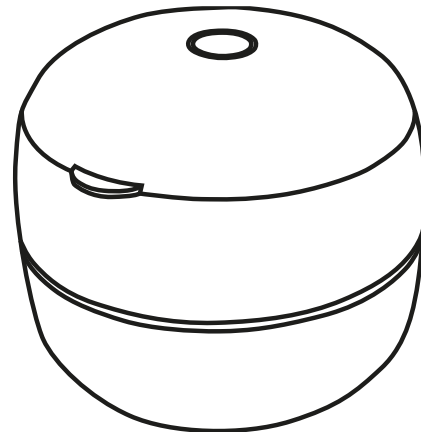


Keep Health Always Around You!

Operation guide

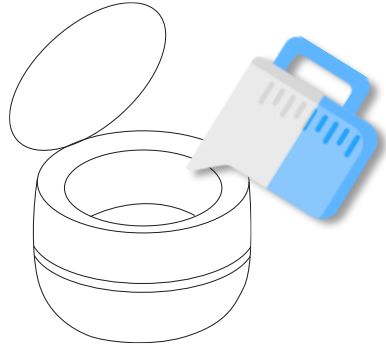


1

Open the lid and fill the stainless steel tank with water.
WARNING: Running a device without water may damage it!

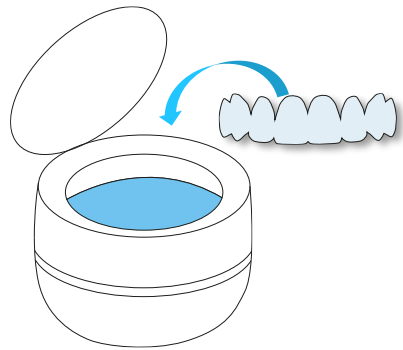

2


Put the objects in the water.
The object should be completely immersed.
Do not exceed the "MAX" level.

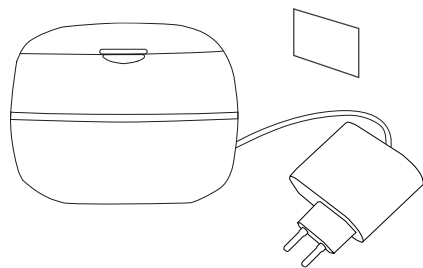

3

Close the lid, and connect the power cord to an outlet.


Operation guide

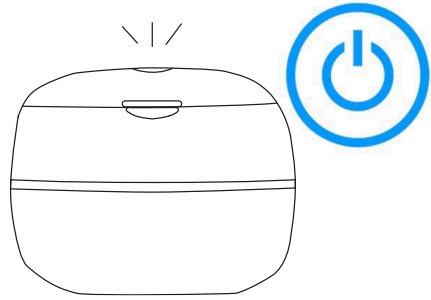

4


Press the button once, and the blue indicator lights up, to start 5 minutes of working time. The indicator will turn off and the cleaner stop working automatically after the duration.

※ If you would like to stop cleaning within 5 minutes of working time, simply press the button again.

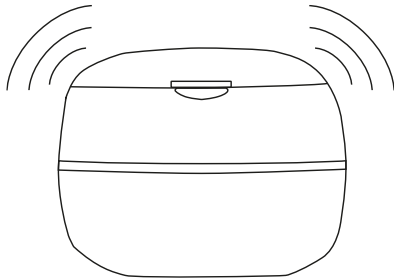

5

While cleaning, a buzzing sound will be heard. This means the device is working.

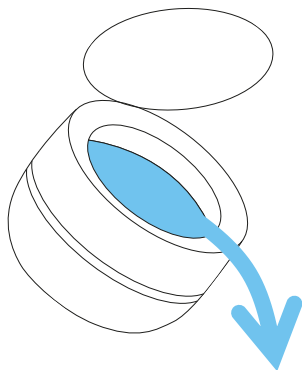



6

After the cleaning process is finished, unplug the power cord, open the lid, and remove the aligner. Lastly, pour the water out of the tank from the arrow direction. Clean and wipe the tank.

※ Note: Do not wet the operation button when pouring out the water.

6-1

When disposing of water, be careful not to let water touch the operation button and power cord.
(Wrap the power cord as shown in the picture to prevent water from entering.)

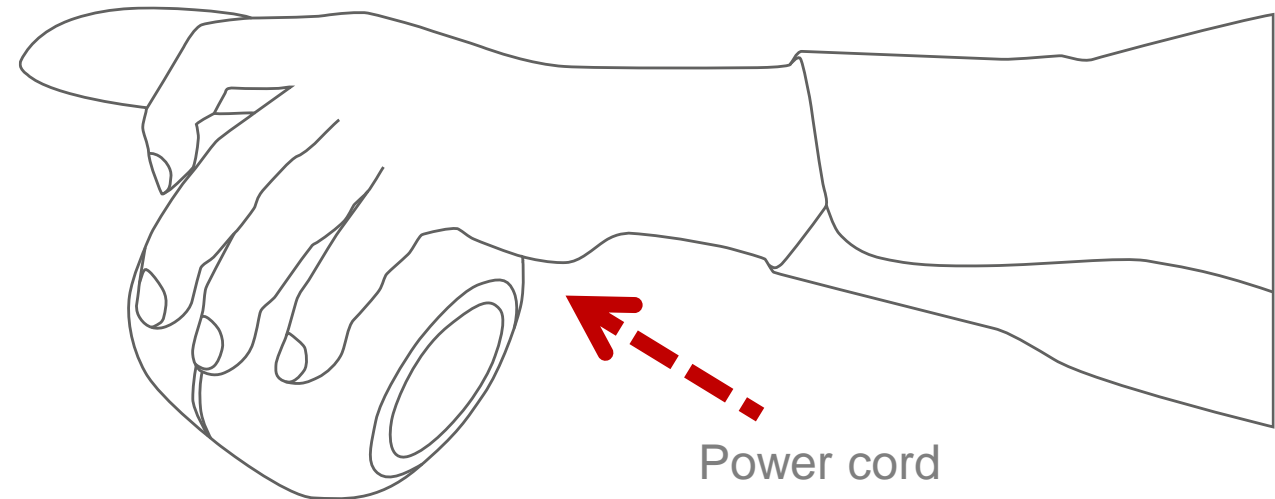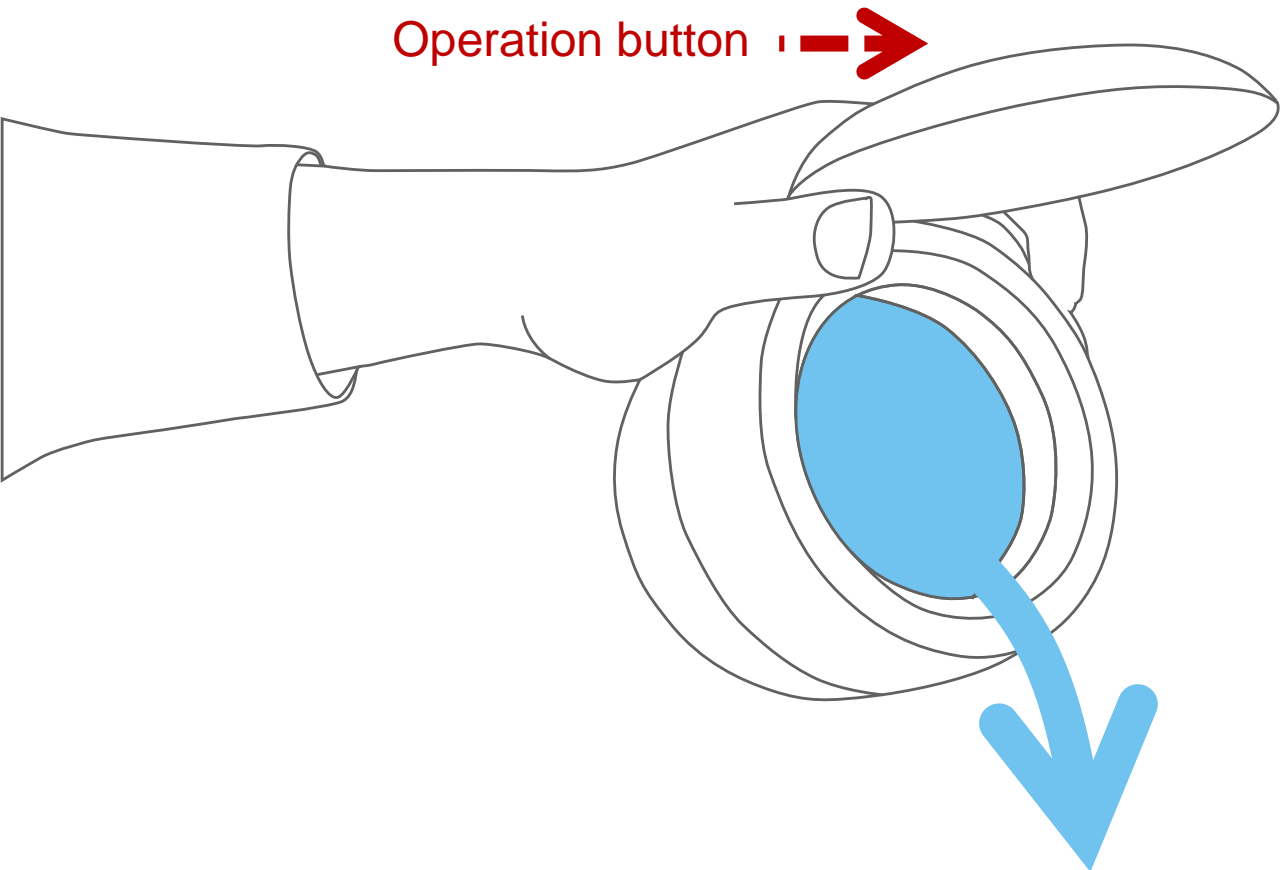



Maintenance

Tank cleaning

The interior needs periodic cleaning.


After use, wipe it with cotton or tissue paper to keep the tank clean and dust-free.

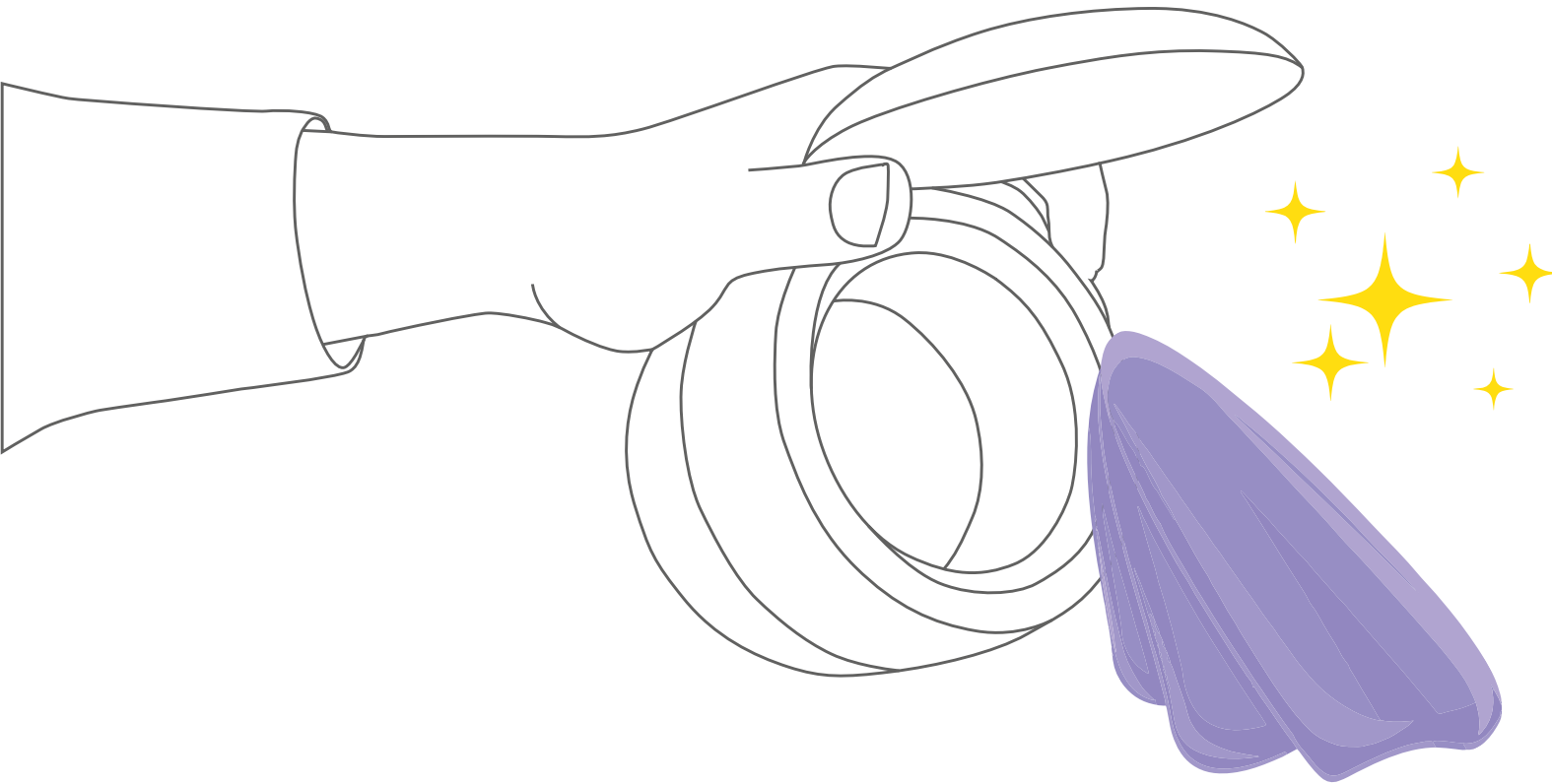

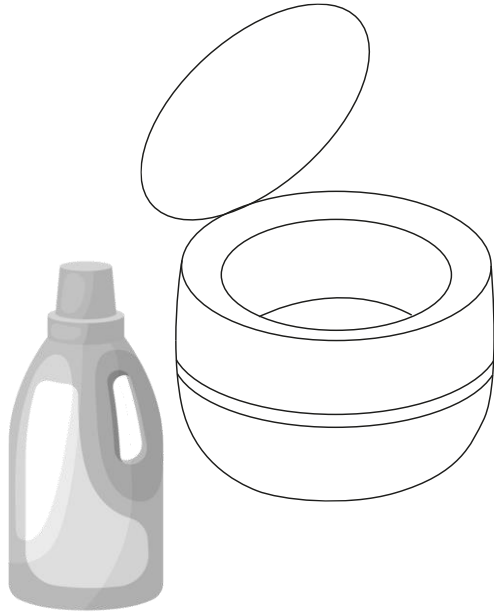

Tank cleaning

- Use water for general cleaning.
 - Thorough cleaning In order to obtain better results, you may add a few drops of universal detergent.
 - We recommend using commercial dishwashing detergent-preferably without lipid-replenishing substances ("balm").
 - For professional cleaning purposes, special cleaning detergents can be used. Find out which concentrate is suitable for your objects to be cleaned and the type of contamination.
- ※ Before cleaning the device, pull out the power plug.
Never immerse the device, cable, or power plug in water or another liquid for cleaning and care purposes.

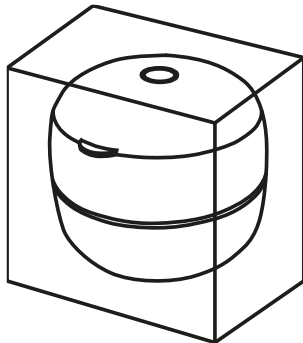

Storage

- Store the device in a cool and dry place.
- Keep the device out of children's reach.



Universal use

In addition to the Aligner, you can clean the items below.

- Personal items: Glasses, Dentures, Toothbrushes, razor blades, baby items, etc.
- Jewelry: Necklaces, rings, earrings, bracelets, watch chains, coins, etc.
- Metal parts: Coins, medals, badges, valves, etc.

※ Loosely framed stones, glued ornaments on decorative items, and loose items could be detached during the cleaning process.

Ask your dentist, particularly in the case of dentures including veneers.



Caution

- When filling the tank, always ensure that the "MAX" marking is not exceeded.
- Do not use aggressive or abrasive cleaners.
- Do not expose the device to temperatures below -10°C or above 40°C.
- CLEANING DETERGENT:

Do not use flammable or explosive detergent.

Special cleaning detergent should be used in case of difficult and/or stubborn dirt, such as grease, resin, rust, wax, pigments, etc. that cannot be removed with water and an all-purpose cleaner.



A DANGER FOR CHILDREN

This device is not intended for use by individuals with restricted physical or mental capacities or those with a lack of experience or knowledge, including children unless they are supervised by an individual responsible for their safety or have received training in operating the device.



- Children shall not carry out cleaning and user maintenance without supervision.
- Store the device and accessories where it is unreachable by children.
- Children must not play with packaging material.



DANGER OF ELECTRIC SHOCK AS A RESULT OF MOISTURE



- Do not use this device in the vicinity of bathtubs, showers, washbasins, or other vessels containing water.
- Do not touch the device with wet hands.
- Prevent the device, the cable, and the mains plug from falling into water or other liquids.
- If the device has fallen into the water, pull out the mains plug immediately, and only then remove the device from the water.



A DANGER OF ELECTRIC SHOCK

- Do not operate the device if it or the cable shows visible signs of damage.
- Only connect the mains plug to a properly installed socket, which is easily accessible and supplies voltage matching the rating plate.
- The socket must be easily accessible after connection.
- Ensures that the cable must not be damaged by sharp edges or hot surfaces.

SPECIFICATION

Name	Tera Harz Care
Info	Ultrasonic clean and warmer for aligners
Overall Size	115 x 115 x103mm, Φ115 x 103H
Tank Size	80 x 80 x 40 mm
Tank Capacity	200ml. Enough space for 2 aligners
Weight	325g
Time	5 min (Automatically turn off)
Max. temp	50°C (After 2 operation,45°C)
Operation Mode	1. Ultrasonic (Once) 2. Warmer + Ultrasonic (Twice) 3. Turn Off (Third)
Ultrasonic Power	24W
Ultrasonic Frequency	42kHz
Input Voltage	AC100-240V
Rated Voltage	12V/2A